

AIRCAST®

ActyFoot™

TAKE CONTROL OF STABILITY



enovis™



AIRCAST®

For more than 30 years, healthcare professionals around the world have trusted Aircast® to provide technological breakthroughs in patient care. Utilizing patented technology supported by biomechanical testing, each Aircast product is developed using the concept of “functional management.” The results are orthopaedic appliances that not only address injuries, but also allow patients to maintain a degree of activity in their daily lives.

Now, following in the footsteps of recognized names such as the AirSelect™ pneumatic walker boot and the Air-Stirrup® ankle brace, Aircast introduces its latest product offering.

ActyFoot™

TAKE CONTROL OF STABILITY.

The ActyFoot™ ankle support's modular design allows the user to control the level of ankle stability as required.

Its removable lateral stay offers increased immobilization and protection to aid recovery following injury or for added confidence during activity.

Without the stay, the unique strap system provides sufficient support for sports and daily activities while allowing greater freedom of movement.



Intended use/Indications

ActyFoot is designed to provide support, stabilization and partial immobilization (restriction of movement) of the ankle

- Acute lateral ligament injuries
- Acute ankle sprains and ligament ruptures
- Prevention of further injuries in pre-operative chronic ankle instability
- Protection of ankle joint
- Functional rehabilitation of ankle injury

Features and Benefits



MODULAR DESIGN

- Detachable lateral stay allows the brace to be modified for use in both acute and chronic injury phases



REMOVABLE SEMI-RIGID LATERAL STAY

- Helps provide additional support and protection by limiting ankle inversion (lateral twisting)



ACCESSIBLE DESIGN

- Can be fully opened to allow simple and secure application



FOREFOOT STRAP

- Allows the wearer to adjust the level of compression as required



CALF STRAP

- Adjustable strap with elastic and non-elastic sections provides additional stability to the ankle joint



CROSS STRAP

- Non-elastic diagonal strap further helps to stabilize the ankle and prevent supination of the foot



SEMI-RIGID FOOT PLATE

- Enhances ankle control without limiting plantar or dorsiflexion



FITTING SYMBOLS

- The symbols on the straps guide the user and facilitate application of the brace



LIGHTWEIGHT, ERGONOMIC

- Anatomically shaped using light materials to help avoid pressure points and provide greater comfort



SLIMLINE DESIGN


- Fits in most traditional athletic shoes or sneakers and can be worn over a sock for additional comfort
- Low profile, thin straps, and durable materials offer a comfortable fit






How it works

Three integrated straps provide proprioceptive feedback and work together with semi-rigid elements to stabilize the ankle and restrict excessive range of movement.

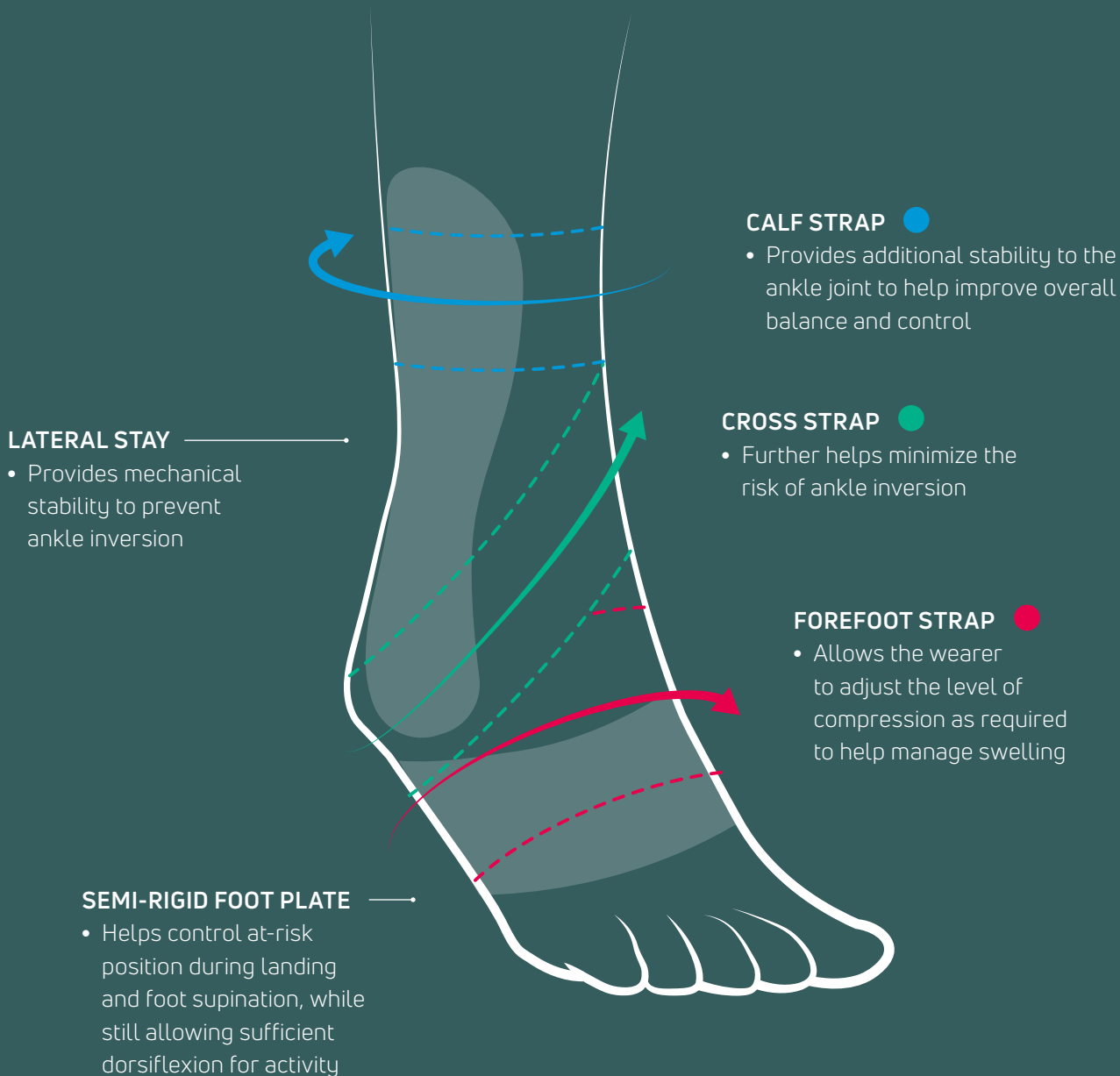

ORDERING INFORMATION

PART NUMBER		HEEL CIRCUMFERENCE*	EU SHOE SIZE	SIZE
LEFT	RIGHT			
82-02A-SL	82-02A-SR	≤ 30 cm	≤ 38	Small
82-02A-ML	82-02A-MR	30 - 34 cm	38-43	Medium
82-02A-LL	82-02A-LR	34 - 36 cm	≥ 43	Large
82-02A-PL	82-02A-PR	≥ 36 cm	≥ 43	Large+

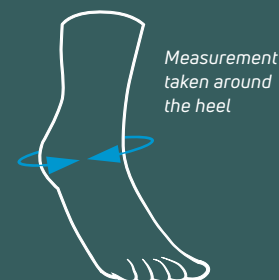

*For more accurate sizing

AIRCAST®

Asia Pacific

1905, Tower 2, Grand Central Plaza
Shatin
HONG KONG
T: 852 31051415
E: marketing.ap@enovis.com

Australia – New Zealand

PO Box 6057
Frenchs Forest DC
NSW 2086
AUSTRALIA
T: 1300 66 77 30
E: customerservice.au@enovis.com

Benelux

Kleinhoefstraat 5, bus 39
Geel 2440
BELGIUM
T: +32 (0)14 248090
E: orders.bx@enovis.com

Canada

6485 Kennedy Rd
Mississauga
Ontario, L5T2W4
CANADA
T: 1-866-866-5031
E: canada.orders@enovis.com

China

A312 SOHO ZhongShan Plaza
1055 W. ZhongShan Rd
Shanghai, 200051
CHINA
T: (8621) 6031 9989
E: shujun.tang@enovis.com

France

3 rue de Bethar
Centre Européen de Frêt
64990 Mouguerre
FRANCE
T: +33 (0)5 59 52 86 90
E: sce.cial@enovis.com

Germany

Bötzingen Straße 90
79111 Freiburg
GERMANY
T: +49 (0) 180 1 676 333
E: kundenservice@enovis.com

Italy

Via Leonardo da Vinci 97
Trezzano sul Naviglio
Milan 20090
ITALY
T: +39 02 835 98001
E: vendite@enovis.com

Nordic

Murmansgatan 126
212 25 Malmö
SWEDEN
T: +46 40 39 40 00
E: info.nordic@enovis.com

South Africa

Unit 6, 2 on London, 2 London Circle
BrackenGate Business Park
Brackenfell, 7560
SOUTH AFRICA
T: +27 (0) 21 276 2968
E: info.southafrica@enovis.com

Spain

Av.Cornella 144. 1-4
08950 Esplugues de Llobregat
Barcelona
SPAIN
T: +34 93 480 32 02
E: ventas@enovis.com

UK & Ireland

1a Guildford Business Park
Guildford, Surrey
GU2 8XG
UK
T: +44 (0)1483 459 659
E: ukorders@enovis.com

USA

2900 Lake Vista Drive
Lewisville
TX 75067
U.S.A.
T: 760 727 1280
E: customercare@enovis.com

enovisTM

Creating Better Together™

enovis-medtech.eu

1a Guildford Business Park, Guildford, Surrey, GU2 8XG, UNITED KINGDOM

Individual results may vary. Neither Enovis, DJO, LLC or any of their subsidiaries dispense medical advice. The contents of this document do not constitute medical, legal, or any other type of professional advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.

MKT-000-00011512-EN - Rev B

Copyright © 2023 DJO, LLC, a subsidiary of Enovis Corporation