

AIRCAST®

ActyToe™

Hinged splint to help bunion pain relief and to aid correction, for use during normal daily activity, during rest and throughout the night.



 **DJO®**

Intended use/Indications

- Pain relief in patients with symptomatic Hallux Valgus (bunion)
- Mild or moderate mal-positioning of the big toe
- Maintaining position after Hallux Valgus (bunion) surgery

Characteristics

- Low-profile design
- Locking clip for splint immobilization (0° and 30° dorsi- and plantar flexion)
- Memory foam cushion
- Adjustable correction angle (0° and 10°)
- Adjustable, thin and soft straps
- Anatomically contoured splint
- Breathable and lightweight materials
- For day and night wear
- Suitable for walking
- Available in 3 sizes (S, M, L)
- For left and right foot
- Can be worn with wide footwear
- Not made with natural rubber latex

Ordering Information

Part Number	Shoe size			Size
	EU	UK Men	UK Women	
82-04S	< 36	< 4	< 3.5	S
82-04M	35 - 42	3.5 - 8	3 - 7.5	M
82-04L	41 - 46	7.5 - 11.5	7 - 11	L

Features and benefits

Available in 3 sizes

For optimal fit

Adjustable mid-foot strap


To support the transverse arch and help stabilize the mid-foot

Articulated splint

Helps achieve pain-free movement of the big toe

Adjustable big toe strap

To help gently correct misalignment of the big toe

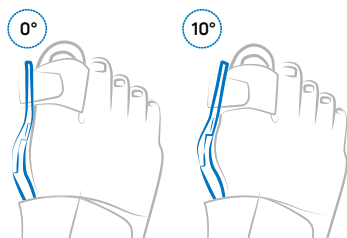


Interchangeable big toe part

- Offers gentle correction at 2 angles (0° or 10°) to aid different degrees of deformation
- For individualized fit and increased patient comfort and compliance

Soft "memory foam" cushion

- To comfort and help protect the bunion
- Aid relief of the metatarsal-phalangeal joint against pressure points

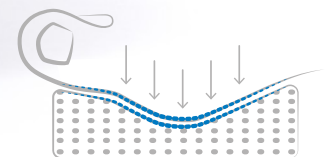

Lockable joint in 0° or 30° dorsi- and plantar flexion

Customized application: three optional locking positions to hold the big toe at a predetermined angle, depending on treatment requirements



DJO Global | 1a Guildford Business Park | Guildford | Surrey | GU2 8XG | UNITED KINGDOM

www.djoglobal.eu

Individual results may vary. Neither DJO Global, Inc. nor any of its subsidiaries dispense medical advice. The contents of this brochure do not constitute medical, legal, or any other type of professional advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.